
Mediblink®

medical devices

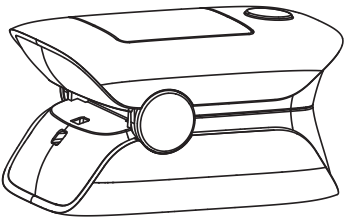EN

Mediblink Fingertip Pulse Oximeter M172

INSTRUCTIONS FOR USE

Please read the instructions for use carefully before using the product.

General Description

Oxygen binds to hemoglobin in red blood cells in the lungs. It is then transported throughout the body in arterial blood. A pulse oximeter uses two wavelengths of light (red and infrared) to determine the percentage of hemoglobin in the blood that is saturated with oxygen. This percentage is called blood oxygen saturation, or SpO₂. A pulse oximeter also measures and displays the pulse rate while measuring the SpO₂ level.

This device conforms to IEC60601-1, IEC60601-1-2, IEC60601-1-11, ISO 80601-2-61.


Measurement Principle

The principle of the pulse oximeter is as follows: The device operates by detecting pulsatile blood flow in the arteriolar vascular bed. The sensor contains a dual light source and a photodetector. One wavelength is 660 nm (red light), and the other is 905 nm (infrared light). Skin, bone, tissue, and venous blood absorb a relatively constant amount of light over time. The photodetector in the sensor collects and converts the transmitted light into an electronic signal proportional to the light intensity. The arteriolar blood flow pulsates and causes variable absorption of light during systole and diastole, as blood volume increases and decreases. The ratio of light absorption during systole and diastole is used to calculate oxygen saturation. This measurement is referred to as SpO₂.

Diagram of Operation Principle

1. Red and infrared-ray emission tube

2. Red and infrared-ray receipt tube

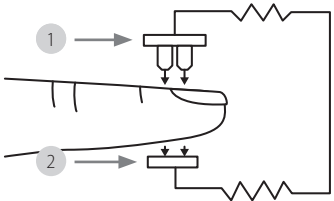Safety Notice

1. Before use, carefully read the manual.

2. Do not use the pulse oximeter:

- if you are allergic to rubber products.
- if the device or finger is damp.
- during MRI or CT scans.
- while taking a blood pressure measurement on the same arm.
- on finger with nail polish, dirt, coatings, or false nail.
- on finger with anatomical changes, edema, scars or burns.
- on a finger that is too large: the finger width is greater than 20mm and the thickness is greater than 15 mm.
- on a finger that is too small: the finger width is less than 10mm and the thickness is less than 5 mm.

• on persons under 18 years of age.

• in conditions with strong changes in environmental light.

• near flammable or explosive gas mixtures.

3. Extended use may cause pain in people with circulatory disorders. Do not use the pulse oximeter for longer than two hours on one finger.

4. Measurements are for informational purposes only. They are not a substitute for a medical examination. If an unexpected reading occurs, the operator may take several additional measurements and consult a doctor.

5. Check the pulse oximeter regularly before use to ensure that there is no visible damage and that the batteries are sufficiently charged. In case of doubt, do not use the device and contact customer service or an authorized retailer.

6. Do not use any additional parts that are not recommended by the manufacturer.

7. Under no circumstances should you open or repair the device yourself. Failure to comply will result in the warranty being void. For repairs, please contact customer service or an authorized retailer.

8. Do not look directly inside the housing during measurement. The red light and the invisible infrared light emitted by the pulse oximeter may be harmful to your eyes.

9. This device is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they are supervised by a person responsible for their safety or have received instructions on how to use the device. Children should be supervised to ensure that they do not play with the device.

10. If the unit has been stored at temperatures below 0 °C, leave it in a warm place for about two hours before use.

11. If the unit has been stored at temperatures above 40 °C, leave it in a cool place for about two hours before use.

12. The displays for the pulse waveform and pulse bar do not allow evaluation of pulse strength or circulation at the measurement site. They are used only to display the current visual signal variation and are not intended for diagnostic purposes.

13. Operation of the fingertip pulse oximeter may be affected by the use of an electrosurgical unit (ESU).

14. Follow local regulations and recycling instructions regarding the disposal or recycling of the device and its components, including batteries.

15. This equipment complies with IEC 60601-1-2 for electromagnetic compatibility for medical electrical equipment and systems. In healthcare facilities or other environments, radio transmission equipment and electromagnetic interference may affect the performance of the oximeter.

16. This equipment is not intended for use during patient transport outside a healthcare facility.

17. When the signal is not stable, the reading may be inaccurate. Please do not rely on the reading.

18. Portable and mobile RF communications equipment can affect medical electrical equipment.

19. WARNING: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

20. WARNING: PORTABLE RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the Fingertip Pulse Oximeter, including cables specified by the MANUFACTURER. Otherwise, degradation of the performance of this equipment may occur.

21. Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Important Testing Guidelines

1. Non-observance of the following instructions may lead to incorrect or failed measurements.

- There must be no nail polish, artificial nails, or other cosmetics on the finger to be measured.
- Ensure that the fingernail is short enough so that the fingertip fully covers the sensor element in the housing.
- Keep your hand, finger, and body steady during the measurement.
- In cases of carbon monoxide poisoning, the pulse oximeter will display a measurement value that is too high.

• To avoid incorrect results, there should be no strong light sources (e.g. fluorescent lamps or direct sunlight) in the immediate vicinity of the pulse oximeter.

• Protect the pulse oximeter from dust, shocks, moisture, and explosive environments.

• Avoid excessive patient movement.

2. The following situations may cause inaccurate measurements:

- Significant levels of dysfunctional hemoglobin (such as carboxyhemoglobin or methemoglobin).
- Venous pulsations.
- Placement of the sensor on an extremity with a blood pressure cuff, arterial catheter, or intravascular line.
- The patient has hypotension, severe vasoconstriction, severe anemia, or hypothermia.
- The patient is in cardiac arrest or in shock.
- Weak pulse quality (low perfusion).
- Low hemoglobin levels.

Intended purpose/Indications/Contraindications

1. Indications for use/Intended Use: The fingertip pulse oximeter is a portable, non-invasive device intended for spot-check measurement of arterial hemoglobin oxygen saturation and pulse rate in adults at home and in hospitals (including clinical use in internal medicine, surgery, anesthesia, etc.).

2. Medical indications: Normal oxygen saturation and hypoxemia.

3. Name of disease or condition: Well-being and diseases that may result in hypoxemia.

4. Patient population: Adults.

5. Intended user: Professionals and lay persons.


6. Contraindications: The device is not for continuous monitoring.

Unit Illustration


Package contents:

- 1 x Pulse oximeter M172 (model XM-111)
- 1 x Owner's Manual
- 1 x Retaining strap
- 1 x Storage bag
- 2 x 1.5 V AAA batteries

Monitor Unit

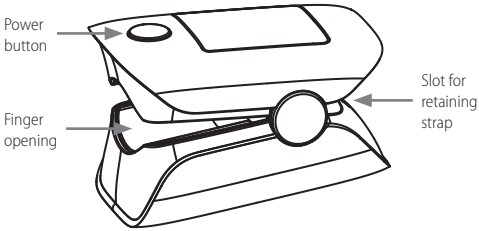Display

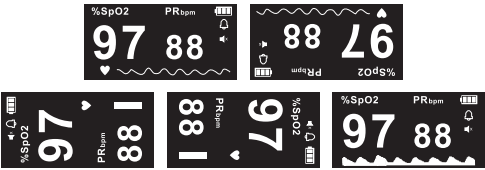Features

1. Simple to operate and convenient to carry.

2. Compact size, lightweight, and low power consumption.

3. Displays SpO₂, PR, pulse bar, and waveform.

4. Adjustable brightness levels (1–5).

5. Five display modes.

6. A low battery warning is displayed when the battery voltage is too low to ensure normal operation of the oximeter.

7. When “Finger out” is displayed, the pulse oximeter will automatically power off after 10 seconds.


8. Audible signal (beep): The device emits a beep synchronized with each detected pulse.

9. When the audible signal is enabled, the displayed values will flash and a distinct audible signal will sound when preset threshold values are exceeded. This audible signal is for informational purposes only and does not constitute an alarm function.

Unit Operation

Battery Installation


Follow the arrow to unlock the battery cover. Install 2 new AAA alkaline batteries according to the correct polarity, install the battery cover and lock it in the opposite direction of the arrow.

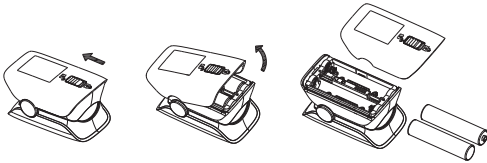Note:

- 1. Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- 2. Use only the size and type of batteries specified.
- 3. Do not mix different types of batteries together or old batteries with new ones. Always replace batteries as a set.
- 4. Replace the batteries in a timely manner when low voltage indicator is lit.
- 5. If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- 6. Do not try to recharge batteries not intended to be recharged; they can overheat and rupture.
- 7. Do not dispose of batteries in fire; batteries may explode or leak.
- 8. Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.
- 9. Please follow local regulations for the disposal of used batteries.

Attaching the Retaining Strap

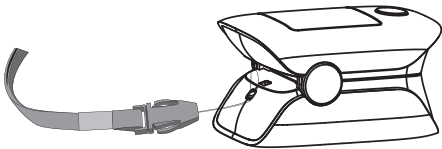System Settings

With the device powered off, press and hold the power button for approximately 2 seconds to enter the settings mode.

Available settings include: Tips, Beep, Language, Default, SpO₂ Tips Lo, PR Tips Hi, PR Tips Lo, and Exit. Press and hold the button to enter a specific setting. Short press the button to switch between setting items.

Settings		Settings	
Tips	* On	SpO2 Tips Lo	* 90
Beep	On	PR Tips Hi	100
Language	EN	PR Tips Lo	50
Default		Exit	

To Use

CAUTION: Ensure that the finger size is appropriate (fingertip width: approximately 10–20 mm; thickness: approximately 5–15 mm).

CAUTION: Do not use this device in a strong electromagnetic radiation environment.

CAUTION: Do not use this device simultaneously with other medical or non-medical devices.

CAUTION: When inserting the finger, ensure it completely covers the LED window inside the finger clip compartment.

1. As shown in the figure, press the clip of the pulse oximeter, fully insert your finger into the finger clip compartment, and then release the clip.

2. Press the power button once on the front panel to turn on the pulse oximeter.

3. Keep your hand still during measurement. Do not move your finger during the test. It is recommended that you remain still while taking a measurement.

4. Read the results on the display screen.

5. To adjust the display brightness, press and hold the power button during operation until the brightness level changes.

Troubleshooting Guide

Problem	Cause	Solution
Monitor does not display.	Batteries are depleted.	Replace the batteries.
	Batteries are not inserted correctly.	Reinsert the batteries. If measurement values are still not displayed after reinserting the batteries correctly, contact customer service.
Measurements are erratic.	Insufficient circulation in the measured finger.	Observe the Important Guidelines.
	Finger, hand or body is moving.	Keep your finger, hand and body still during the measurement.
	Cardiac arrhythmia.	Seek medical attention.
Measurements cannot be displayed normally.	Finger is not inserted correctly.	Retry by inserting the finger correctly.
	Patient's SpO ₂ value is too low to be measured.	There may be excessive ambient light. Try again several times. If you are sure there is no problem with the product, go to a hospital promptly for an accurate diagnosis.

Specification

Model/REF	XM-111/M172	
Product	Fingertip pulse oximeter	
Display	0.96" display	
	Display Range	0%–100%
SpO ₂	Measurement Range	70%–100%
	Accuracy	70%–100% ±2%
		0%–69% not defined
	Resolution	1%

Pulse Rate	Display Range	0–240bpm
	Measurement Range	30–240bpm
	Accuracy	30–100bpm, ±2bpm; 101–240bpm, ±2%
	Resolution	1bpm
Power supply	2 × 1.5 V AAA batteries	
Weight	Approx. 54g	
Dimensions	Approx. 60.2 × 35 × 35.5 mm	
Operating Environment	Temperature	5°C–40°C
	Humidity	15%–93%RH
	Pressure	700hPa–1060hPa
Storage Environment	Temperature	–20°C–55°C
	Humidity	15%–93%RH
	Pressure	700hPa–1060hPa
Ingress Protection Rating	IP22	
Classification	Internal Powered Equipment Type BF	
Applied parts	Finger pad	
The Data UPDATE period		Less than 12s

Note: The functional tester cannot be used to assess the accuracy of the oximeter. The methods used to determine SpO₂ accuracy are based on clinical testing. The oximeter measures arterial haemoglobin oxygen saturation (SpO₂), and the results are compared with values obtained from arterial blood samples analysed using a CO-oximeter.

- EN ISO 80601-2-61: Medical electrical equipment – Part 2-61: Particular requirements for the basic safety and essential performance of pulse oximeter equipment.
- EN 60601-1: Medical electrical equipment – Part 1: General requirements for basic safety and essential performance.
- EN 60601-1-2: Medical electrical equipment – Part 1-2: General requirements for basic safety and essential performance – Collateral standard: Electromagnetic disturbances – Requirements and tests.
- EN 60601-1-11: Medical electrical equipment – Part 1-11: General requirements for basic safety and essential performance – Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment.

Correct disposal of this product
(Waste Electrical and Electronic Equipment – WEEE)

This symbol on the product indicates that it must not be disposed of with other household waste at the end of its service life. To prevent potential harm to the environment or human health, this product should be disposed of separately and recycled responsibly. For household users: please contact the retailer from whom the product was purchased or your local government office for details on environmentally safe recycling options. For business users: please contact your supplier and check the terms and conditions of the purchase agreement. This product must not be mixed with other commercial waste for disposal. This product does not contain hazardous substances.

Icon Explanation

Symbol	Definition	Symbol	Definition
	CE Marking of Conformity		Humidity limitation
	Conformity to WEEE Directive		Atmospheric pressure limitation
	Caution		Serial number
IP22	Ingress protection degree		Alarm inhibit
	Direct Current	PR _{bpm}	Pulse rate (BPM)
	Type BF applied part		Medical device
	Refer to instruction manual		Unique device Identifier
http://www.mediblink.com/f/m172.pdf			Non-ionizing electromagnetic radiation
	Authorized Representative in European Union		Low power indication
SpO₂%	Oxygen saturation		Beep OFF
	Manufacturer		Beep ON

Symbol	Definition	Symbol	Definition
	Date of Manufacture		Tips OFF
	Batch Number*		Tips ON
	Temperature limit		Number of products included in packaging
	Paper		

* The production date can be read from LOT number [YYYYMM]; first four digits represent the Year and last two digits Month of production. Example: LOT 202503 = March 2025

Statement of Essential performance

- When the oximeter is placed on the patient's finger or a simulated finger, SpO₂ and PR values are displayed normally.
- Measurement accuracy:
 - Clinical accuracy of SpO₂: within the range of 70%–100% ($Arms \leq 3\%$);
 - For each specified range, the SpO₂ accuracy of the pulse oximeter shall be expressed as the root mean square (*rms*) difference between the measured values (SpO_{2i}) and the reference values (*S_{Ref}*).

$$Arms = \sqrt{\frac{\sum_{i=1}^n (SpO_{2i} - S_{Ref})^2}{n}}$$

Probe Accuracy of Oximeter

Items	Descriptions
Clinical SpO ₂ accuracy (70–80%)	Arms: < 3
Clinical SpO ₂ accuracy (80–90%)	
Clinical SpO ₂ accuracy (90–100%)	
Sterile	Non-sterile
Default settings	None
Reuse	Disinfect for repeated use

- Instructions for inspection of the application site (skin integrity): Before inserting a finger into the oximeter probe, visually inspect the skin to ensure it is intact and free from injury or other abnormalities.
- Instructions regarding sensor relocation: No sensor relocation is required during the service life of the product.
- Use during movement and low perfusion conditions:
 - Do not move your finger, arm, or body during measurement. Movement, including talking, coughing, or sneezing, may affect the accuracy of the measurement results.
 - Measurements may be unreliable under conditions of low perfusion.

The sensor and monitor have been validated and tested for compliance with EN ISO 80601-2-61:2017.

Clinical data were obtained from a controlled induced hypoxia study conducted on healthy adult volunteers. A total of 12 subjects (5 females and 7 males), aged 22 to 45 years and non-smokers, were included. In the clinical evaluation, 288 paired data points from the 12 subjects were analysed.

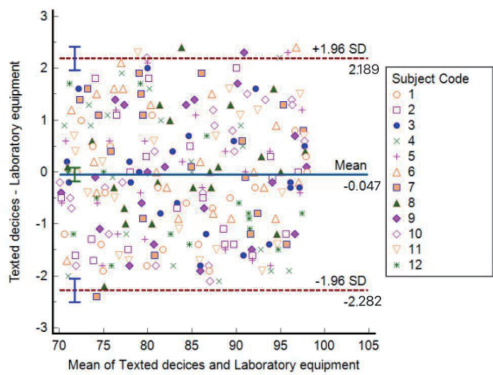
Fig.1: Bland-Altman plot of SaO₂ versus (SpO₂-SaO₂)

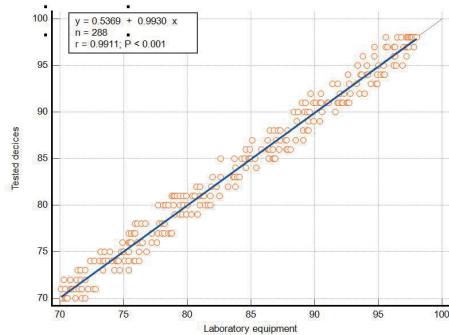
Fig. 2: Regression plots for all subjects

Electromagnetic Compatibility Information

The device complies with the EMC requirements of the international standard IEC 60601-1-2. These requirements are met under the conditions described in the table below. This device is medical electrical equipment and requires special precautions regarding electromagnetic compatibility (EMC), as described in these instructions for use. Portable and mobile RF communications equipment may affect the device. The use of non-approved accessories may negatively affect the device and alter its electromagnetic compatibility. The device should not be used directly adjacent to or stacked with other electrical equipment. If such use is necessary, the device should be observed to verify normal operation.

Table 1

Electromagnetic Compatibility Information – electromagnetic emission	
This device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.	
Emissions test	Compliance
CISPR 11 Conducted emissions	Not applicable
CISPR 11 Radiated emissions	Group 1 Class B
Harmonic emissions IEC 61000-3-2	Not applicable
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable

Table 2

Guidance and manufacturer's declaration – electromagnetic immunity		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Immunity test	IEC 60601 test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8kV contact ±2kV, ±4kV, ±8kV, ±15kV air	±8kV contact ±2kV, ±4kV, ±8kV, 15kV air
Electrostatic transient / burst IEC 61000-4-4	±2kV for power supply lines 100kHz repetition frequency ±1kV for input/output lines	N/A
Surge IEC 61000-4-5	±0.5kV, ±1kV differential mode line-line	N/A
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% U _i (100% dip in U _i) for 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°	N/A
	0% U _i (100% dip in U _i) for 1 cycle at 0°	
	70% U _i (30% dip in U _i) for 25/30 cycles at 0°	
	0% U _i (100% dip in U _i) for 250/300 cycle at 0°	
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30A/m, 50/60Hz	30A/m, 50/60Hz

NOTE: U_i is the a. c. mains voltage prior to application of the test level.

Table 3

Guidance and manufacturer's declaration – electromagnetic immunity		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Immunity test	IEC 60601 test level	Compliance level
Conducted RF IEC 61000-4-6	3Vrms 150kHz to 80MHz 6Vrms 150kHz to 80MHz outside ISM bandsa	N/A
Conducted RF IEC 61000-4-3	10V/m 80MHz to 2.7GHz	10V/m

Proximity magnetic fields IEC 61000-4-39	30kHz, 8A/m, CW 134,2kHz, 65A/m Pulse modulation 2,1kHz 13,56kHz, 7,5A/m Pulse modulation 50kHz	30kHz, 8A/m, CW 134,2kHz, 65A/m Pulse modulation 2,1 kHz 13,56 kHz, 7,5A/m Pulse modulation 50kHz
	See Table 5	See Table 5
Proximity fields from RF wireless communications equipment IEC 61000-4-3	See Table 5	See Table 5

Table 4

Recommended separation distances between RF wireless communications equipment			
The device is intended for use in the electromagnetic environment in which radiant RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between RF wireless communications equipment and the device as recommended below, according to the maximum output power of the communication equipment.			
Frequency MHz	Maximum power W	IEC 60601 level	Compliance level
385	1.8	27	27
450	2	28	28
710	0.2	9	9
745			
780			
810	2	28	28
870			
930			
1720			
1845	2	28	28
1970	2	28	28
2450			
5240			
5500			
5785	0.2	9	9

Note 1: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

WARNINGS!

- This device should not be used in the vicinity of, or stacked with, other electrical equipment such as mobile phones, transceivers, or radio-controlled devices. If such use is necessary, the device should be observed to verify normal operation.
- The use of accessories and power cords other than those specified, with the exception of cables sold by the manufacturer as replacement parts for internal components, may result in increased emissions or decreased immunity of the equipment or system.
- The operator shall ensure that the device maintains basic safety and essential performance with regard to EMC.

Table 5

Guidance and declaration of manufacturer – electromagnetic immunity				
Test frequency (MHz)	Band ^{a)} (MHz)	Service ^{a)}	Modulation	IMMUNITY TEST LEVEL (V/M)
385	380 to 390	TETRA 400	Pulse modulation ^{b)} 18Hz	27
450	430 to 470	GMRS 460, FRS 460	FM ^{c)} ±5 kHz deviation 1 kHz sine	28
710	704 to 470	LTE Band 13,17	Pulse modulation ^{b)} 217Hz	9
745				
780				
810	800 to 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation ^{b)} 18Hz	28
870				
930				
1720				
1845	1700 to 1990	GSM 1800, CDMA 1900; GSM 1900; DECT; LTE Band 1,3,4,25; UMTS	Pulse modulation ^{b)} 217Hz	28
1970	2400 to 2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation ^{b)} 217Hz	28
2450				
5240	5100 to 5800	WLAN 802.11 a/n	Pulse modulation ^{b)} 217Hz	9
5500				
5785				

If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1m. The 1m test distance is permitted by IEC 61000-4-3.

- For some services, only the uplink frequencies are included.
- The carrier shall be modulated using a 50% duty cycle square wave signal.
- As an alternative to FM modulation, the carrier may be pulse modulated using a 50% duty cycle square wave signal at 18Hz. While it does not represent actual modulation, it would be worst case.

Warranty

Product: Mediblink Fingertip Pulse Oximeter M172
Manufactured for (Importer for EU & Distributor):
Mediblink d.o.o., Gubčeva cesta 19, 8210 Trebnje, Slovenija,
info@mediblink.com, www.mediblink.com

Seller's name, address, signature and stamp*:

Date of delivery/sales*:

* If the invoice is accompanied by this warranty, and if all above information can be seen from the invoice, it is not necessary to fill in this field.

WARRANTY TERMS

Dear customers!

The warranty period for the device is **2 years** and begins on the date of purchase or delivery of the product to the customer. The warranty does not apply to batteries. The warranty applies to original products purchased from sellers or distributors who provide valid proof of purchase and product traceability. When claiming the warranty, the customer must provide the invoice or other proof of purchase, or a completed warranty certificate. The consumer's statutory rights arising from the conformity of goods remain unaffected.

Many complaints are the result of improper use of the product. You can easily avoid any problem, by getting useful information provided by our special service department. To reach our service department, you can call or send e-mail to Mediblink local distributor.

Before sending the product back to retailer, we kindly ask you to call our service department, to get help about how to use the device to save you with unneeded trips.

The manufacturer guarantees free elimination of all imperfections due to defects in material or manufacturing procedure by repairing or replacing the product. In case that the product can not be repaired or replaced, the customer will get the money refund. The guarantee is not valid in case of force majeure, accidents or unexpected events (such as lightning, water, fire etc.), incorrect use or incorrect transport, non-compliance with safety and maintaining regulations or in case of unprofessional product intervention.

Traces of every day product usage (scratches, abrasions) and not subject to claim. The warranty does not eliminate the customer rights, which originate from seller responsibility for product flaws. By accepting the claimed product by the service department, the service department does not take responsibility for loss of saved data or settings on the product. All product repairs, which are performed out of product warranty period, have to be paid by customer by prior notice.

The manufacturer guarantees the product quality and flawless product operation in the warranty period, which starts with the day of product purchase. If the product can not be repaired in 45 days, the product will be replaced with a new one. In case that the product can not be replaced, the money will be refunded to the customer.

In case of product claim, call or send e-mail to Mediblink local distributor.

Manufactured for (Importer for EU & Distributor):
Mediblink d.o.o., Gubčeva cesta 19,
8210 Trebnje, Slovenija,
info@mediblink.com, www.mediblink.com

Manufacturer:
JOYTECH Healthcare Co., Ltd.
No.365, Wuzhou Road, 311100 Hangzhou,
Zhejiang Province, People's Republic of China
E-mail: info@sejoy.com
 Shanghai International Holding Corp. GmbH (Europe)
Eiffestrasse 80, 20537 Hamburg, Germany
E-mail: shholding@hotmail.com

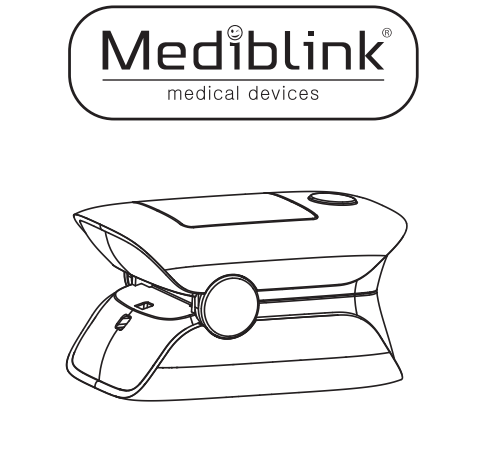
Instructions for use version: E
Date of issue: 02. 09. 2025
Date of last change: 20. 4. 2026

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